


MANUFACTURER AND EXPORTER OF QUALITY
Surgical, Dental, Veterinary Instruments.

All kind of Single Use Surgical instruments
Available in our regular instruments catalogue

edition 2020/21 en



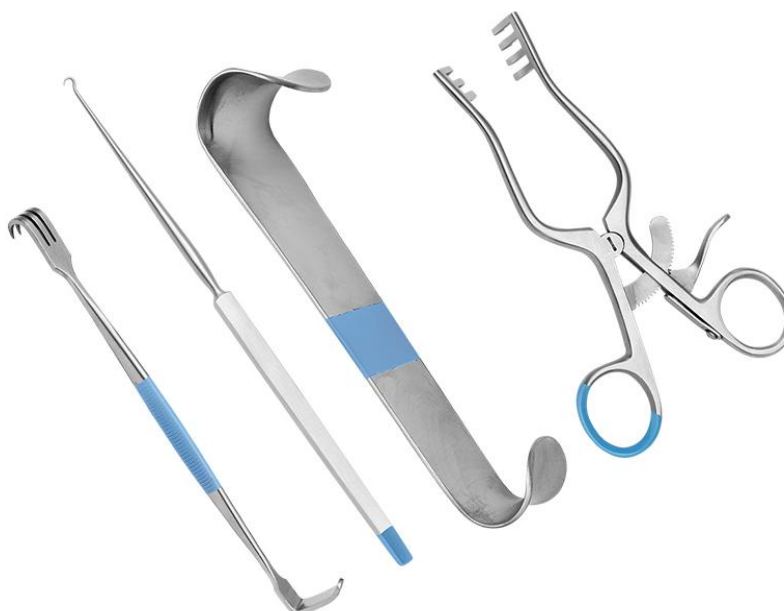


② Single Use surgical instruments




Single use surgical instruments

Single use surgical instruments protect patients From infection. They are also cost effective, Especially when compared to the whole cost of using reusable instruments including initial Purchase, the cost of replacing lost instruments and repairing broken ones and transportation to and from sterilizing centers or the cost of running a desk top sterilizer.



Item Code	FINISH	AVAILABLE COLOURS					NOTES:-	DRAWING#:	Approved by: Q.A Manager	Scale: 1:1	MEAS.UNIT: MM
03-576	Satin						1: All Dimensions are in mm.	Rev	Description	Date	DEVICE:
		<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input type="checkbox"/>	2: Material AISI 420 Stainless Steel	00	Original	16-12-2017	Iris Scissor CVD 11.5 cm
							3: All colors are available in all colors				DESCRIPTION:
											Iris Scissor CVD 11.5 cm
HARDNESS: 43 to 53										16-12-2017	GOLD SKY INT.

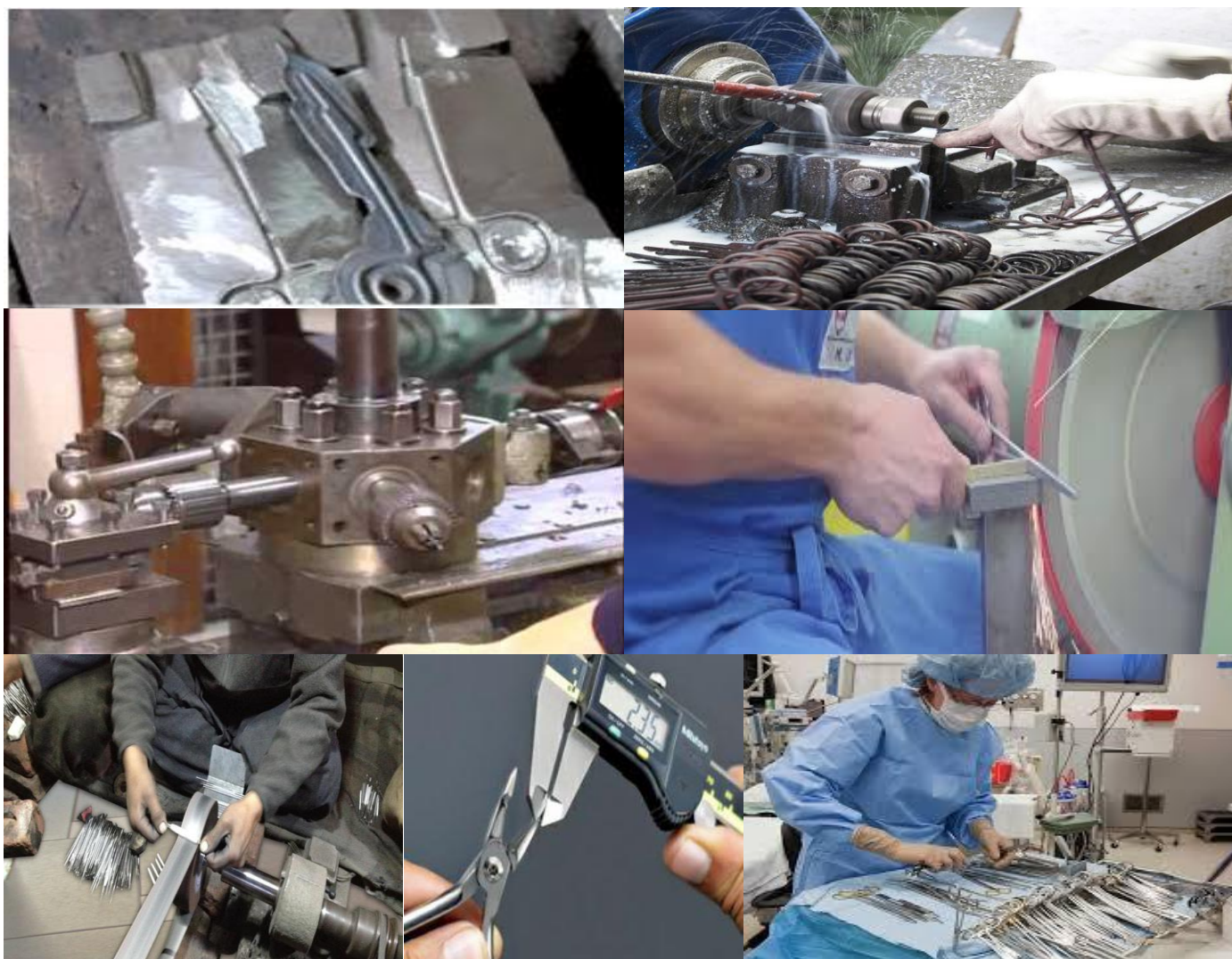

WIDE RANGE OF HIGH QUALITY SINGLE-USE SURGICAL INSTRUMENTS

Thanks to continuous innovation and the use of certified raw materials, GOLD SKY INTERNATIONAL has succeeded in creating a complete range of instruments Capable of satisfying the most sophisticated needs of an international audience both in terms of quality and range.

Thanks to good relationship with clients and satisfaction of our clients. To Date GOLD SKY INTERNATIONAL is present with its own branch in Italy.

GOLD SKY INTERNATIONAL, IS BASED IN SIALKOT, WHICH IS FAMOUS FOR ITS INDUSTRY OF, SURGICAL, SPORTS, LEATHER. WE WELCOME YOU TO VISIT US OR CONTACT US DIRECTLY. WE SHALL PROVIDE YOU OUR BEST SERVICES.

MANUFACTURING PROCESS

FORGING:

- .SELECTION OF FORGING PROCESS. i.e. , hot ,warm and cold forging.
- .Forging process apply on instruments/devices.(i.e. hot, warm and cold) for shaping metal parts.
- .Trimming of excess Material and acid cleaning.

MILING AND TURNING

- .Cutting away unnecessary material by milling machines.
- .Making serrations, Ratchets, male and female box joint and grooves.

FILING, GRINDING AND SETTING

- .Removing of Irregularity by grinding wheels and files so instruments develop its final shap.
- .Initial setting of instruments by hand with different type of gauges.

HEAT TREATMENT

- .Heat treatment is significant to improve performance and strength of instruments.
- .Heat treatment process is done to change physical and mechanical properties of instruments.
- .Without modification of original shape of instruments
- .For soften the material / metal annealing process done
- .Hardening and tempering process to done to produce durability in instruments.
- .Rock well hardness (HRC) for needle holders without tungsten carbide inserts around 40-48 while Scissor among HRC 50-48

ELECTRO-POLISHING

- .The Electro-polishing departments are also managed by trained and skilled workmen using the imported chemicals and grinding wheels to ensure the best- Quality finish.

TUNGSTEN CARBIDE

Our T.C Department is installed with best Quality imported Machinery being managed by the exceptionally well trained Expert using the German Technology.

PASSIVATION

Using the top Quality imported chemical in our passivation department. We passivate and do boil test on each and every instrument before packing – subjecting every instruments 24 hour treatment and if need to be re-passivation process so as to eliminate any chances of rusting/ Corrosion etc.

SAND BLASTING

For sand blasting our instruments. We use the machine as well as the best quality imported sand blasting material which are using by the well trained and skilled operators.

FINAL INSPECTION

To make sure Quality of instruments per desired specifications final checking/inspection done by qualified Technician.

MARKING /LABELING /PACKAGING

- .Marking/printing process done by laser marking system. Etching and screen printing
- . Packaging done per desired specification



Adson tissue forceps 1x2 12cm
GS-5001-12



Adson forceps serrated 12cm
GS-5002-12



Adson Brown forceps 12cm
GS-5003-12 , 7x7



Martin Splinter 11.5cm
GS-5004-11.5



Cushing Dissecting forceps
20cm
GS-5005-20



Dissecting Forceps medium
medium 14cm
GS-5006-14



Dissecting forceps
standard 14cm
GS-5007-14



Debakey Atraumatic
Dissecting Forceps
GS-5008-15



**Lane Tissue forceps
1x2, 1615cm
GS-5015-15**



T.O.E forceps
12.5cm
GS-5009-12.5



Treves forceps
toothed 12.5cm
GS-5010-12.5



Gillies forceps
toothed 15cm
GS-5011-15



McIndoes
non tooth 15cm
GS-5012-15



Iris forceps
toothed 10cm
GS-5013-10



Iris forceps
non tooth 10cm
GS-5014-10



Lister Bandage scissor
14cm
GS-5501-14



Lister Bandage scissor
18cm, 20cm
GS-5502-18



Braun Staddler Scissor
14cm
GS-5503-14



Braun Staddler Scissor
18cm
GS-5504-18



Operating Scissor S/B
14cm
GS-5505-14



Operating Scissor B/B
14cm
GS-5506-14



Operating Scissor S/S
14cm
GS-5507-11.5



Delicate Scissor 9cm
GS-5508-9



Metzenbaum Scissor
14cm
GS-5509-14



Iris Scissor 11.5cm
GS-5510-11.5



Noyes eye Scissor
12.5cm
GS-5511.12.5



**umbilical cord
Scissor 10,5cm
GS-5510**



**Bush Scissor
5518-14
5519-16**



Mayo Hegar N/H 14cm

GS-6001-14



Mayo Hegar N/H 16cm

GS-6002-16



Mayo Hegar N/H 18cm

GS-6003-18



Mayo Hegar N/H 20cm

GS-6004-20



Derf needle holder 12cm

GS-6005-12.5



Halsey needle holder 12cm

GS-6006-12



Mosquito forceps 12.5cm

GS-6007-12



Kocher forceps 14cm
GS-5008-14



Kocher forceps 16cm
GS-5009-16



Kocher forceps 18cm
GS-5010-18



Kocher forceps 20cm
GS-5011-20



Rankin Crile Artery forceps
14cm
GS-5012-14



Rankin Crile Artery forceps
16cm
GS-5014-14



Kelly Artery forceps
14cm str
GS-5016-14



Kelly Artery forceps
14cm cvd
GS-5017-14



Spencer Wells Artery forceps.
13cm
GS-5018-13



Spencer Wells Artery forceps.
15cm
GS-5020-16



Spencer Wells Artery forceps.
18cm
GS-5022-18



Spencer Wells Artery forceps.
20cm
GS-5024-20



Wietlaner Retractor
13cm 3x4 prong
GS-6001-13



Roux Retractor 17cm
18 x 5.5mm
GS-6002-17



Senn Miller Retractor
16.5cm
GS-6003-16.5



Williger Bone curette
GS-6004



**Scalple Handle no.3
GS-6021**



**Bone Curette
double ended
GS-6005**



**Bone Curette
Double ended
16cm
GS-6006**



**Single Use Dissecting kit
set of 6pcs.
GS-7000**



**Allis Tissue forceps
16cm
GS-5026-16**



**Little wood forceps
21cm
GS-5027.21**



**Babcock forceps
16cm
GS-5028**



**Crile wood needle
holder 15cm
GS-6029**



Luer vulsellum 1x1
24cm
GS-5030-24



Allis Tissue Forceps
11"
GS-5031-11



Teals Vulsellum cvd
3x4 teeth, 23cm
GS-5032



Remplay Sponge
24cm
GS-5033-24



Bony polypus
Sponge 23cm
GS-5035-23



Magil forceps Adult
25cm
GS-5036-25



Sims urine curette
Double ended
GS-6008



SINGLE USE INSTRUMENTS



Winterton Speculum
GS-6009



Sims Speculum
GS-6015-S
GS-6016-M
GS-6017-L




CUSCO Speculum
GS-6010-virgin
GS-6011-small
GS-6012-Medium
GS-6013-Large
GS-6014-Medium long




Graves Speculum
GS-6018-S
GS-6019-M
GS-6020-L



**Uerine Hegar dilator
set of 8
GS-6022**

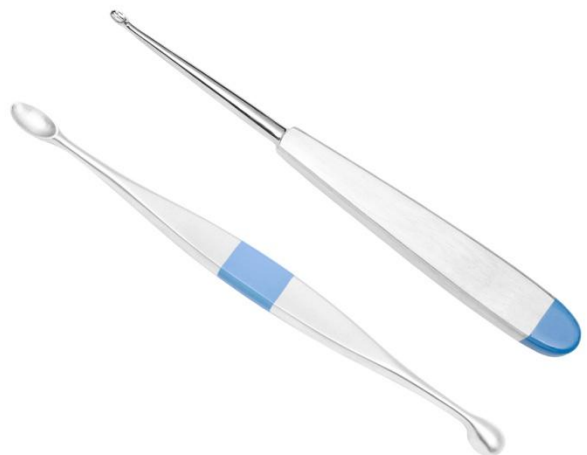

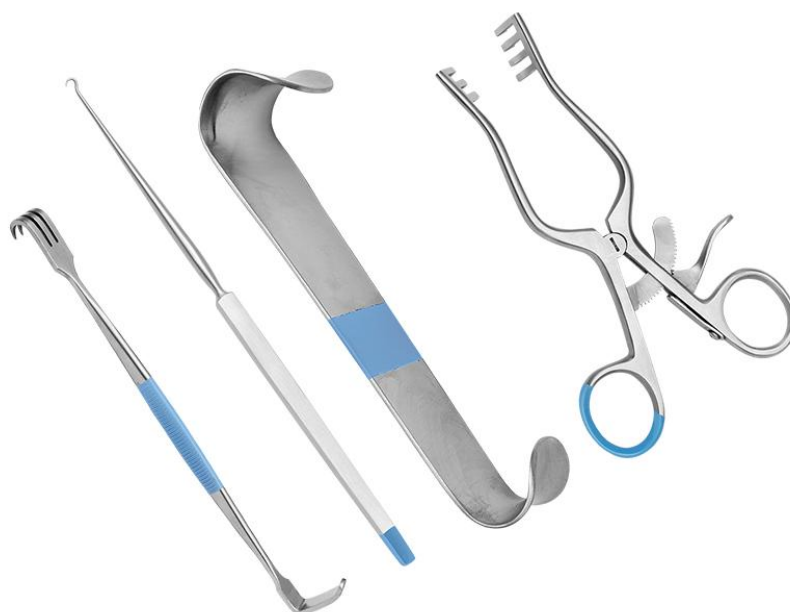



**ISO 9001-20015
ISO 13485-2016
CERTIFIED**



**CE
Approved**

Head office



GOLD SKY INTERNATIONAL

Marala Road, Wakeel Town, Murad Pur.

SIALKOT. PAKISTAN.
TEL: 0092.52.4564095

FAX: 0092.52.4564195
info@goldskyint.com

EU REPRESENTATIVE



Pak Star CO.

VIA REPUBBLICA NO.4,
25069 VILLA CARCINA BRESCIA ITALY.
TELFAX NO: 00390308031585.
info@pakstarco.com
www.pakstarco.com